

### SECTION 1: Identification of the substance/mixture and of the company/undertaking

#### 1.1. Product identifier

Product form : Mixture  
Product name : ZirClean

#### 1.2. Relevant identified uses of the substance or mixture and uses advised against

##### 1.2.1. Relevant identified uses

Use of the substance/mixture : For Rx Only

##### 1.2.2. Uses advised against

No additional information available

#### 1.3. Details of the supplier of the safety data sheet

##### Manufacturer

BISCO, Inc.  
1100 W. Irving Park Rd.  
Schaumburg, IL 60193  
U.S.A.  
T 1-800-247-3368 or 1-847-534-6000  
[www.bisco.com](http://www.bisco.com)

##### EC REP

BISCO France  
208, allée de la Coudoulette  
13680 Lançon de Provence  
France  
T 33-4-90-42-92-92

#### 1.4. Emergency telephone number

Emergency number : CHEMTREC - 24-Hour Hazmat Emergency Communications Center  
Domestic: 1-800-424-9300 Outside the U.S.: 1-703-527-3887, collect calls accepted

### SECTION 2: Hazards identification

#### 2.1. Classification of the substance or mixture

##### Classification according to Regulation (EC) No. 1272/2008 [CLP]

Skin corrosion/irritation, Category 1 H314  
Full text of H- and EUH-statements: see section 16

##### Adverse physicochemical, human health and environmental effects

Causes severe skin burns and eye damage.

#### 2.2. Label elements

##### Labelling according to Regulation (EC) No. 1272/2008 [CLP]

Hazard pictograms (CLP) :



GHS05

Signal word (CLP) :

Danger

Hazard statements (CLP) :

H314 - Causes severe skin burns and eye damage.

Precautionary statements (CLP) :

P260 - Do not breathe fume/mist/vapours.

P264 - Wash hands thoroughly after handling.

P280 - Wear protective gloves/protective clothing/eye protection.

P301+P330+P331 - IF SWALLOWED: Rinse mouth. Do NOT induce vomiting.

P303+P361+P353 - IF ON SKIN (or hair): Take off immediately all contaminated clothing.

Rinse skin with water or shower.

P304+P340 - IF INHALED: Remove person to fresh air and keep comfortable for breathing.

P305+P351+P338 - IF IN EYES: Rinse cautiously with water for several minutes. Remove

contact lenses, if present and easy to do. Continue rinsing.

P310 - Immediately call a POISON CENTER, a doctor.

# ZirClean

## Safety Data Sheet

according to Regulation (EC) No. 1907/2006 (REACH) with its amendment Regulation (EU) 2020/878

P321 - Specific treatment (see supplemental first aid instruction on this label).  
P501 - Dispose of contents and container to a hazardous or special waste collection point, a licensed hazardous-waste disposal contractor or collection site except for empty clean containers which can be disposed of as non-hazardous waste, hazardous or special waste collection point, in accordance with local, regional, national and/or international regulation.

### 2.3. Other hazards

Contains no PBT/vPvB substances  $\geq 0.1\%$  assessed in accordance with REACH Annex XIII

Component	
Potassium Hydroxide (1310-58-3)	This substance/mixture does not meet the PBT criteria of REACH regulation, annex XIII This substance/mixture does not meet the vPvB criteria of REACH regulation, annex XIII

The mixture does not contain substance(s) included in the list established in accordance with Article 59(1) of REACH for having endocrine disrupting properties, or is not identified as having endocrine disrupting properties in accordance with the criteria set out in Commission Delegated Regulation (EU) 2017/2100 or Commission Regulation (EU) 2018/605 at a concentration equal to or greater than 0,1 %

## SECTION 3: Composition/information on ingredients

### 3.1. Substances

Not applicable

### 3.2. Mixtures

Name	Product identifier	%	Classification according to Regulation (EC) No. 1272/2008 [CLP]
Potassium Hydroxide	CAS-No.: 1310-58-3 EC-No.: 215-181-3 EC Index-No.: 019-002-00-8	< 10	Acute Tox. 4 (Oral), H302 Skin Corr. 1A, H314

### Specific concentration limits:

Name	Product identifier	Specific concentration limits
Potassium Hydroxide	CAS-No.: 1310-58-3 EC-No.: 215-181-3 EC Index-No.: 019-002-00-8	( 0.5 $\leq$ C < 2) Eye Irrit. 2, H319 ( 0.5 $\leq$ C < 2) Skin Irrit. 2, H315 ( 2 $\leq$ C < 5) Skin Corr. 1B, H314 ( 5 $\leq$ C < 100) Skin Corr. 1A, H314

Full text of H- and EUH-statements: see section 16

## SECTION 4: First aid measures

### 4.1. Description of first aid measures

First-aid measures general : Call a physician immediately.  
First-aid measures after inhalation : Remove person to fresh air and keep comfortable for breathing.  
First-aid measures after skin contact : Rinse skin with water/shower. Take off immediately all contaminated clothing. Call a physician immediately.  
First-aid measures after eye contact : Rinse cautiously with water for several minutes. Remove contact lenses, if present and easy to do. Continue rinsing. Call a physician immediately.  
First-aid measures after ingestion : Rinse mouth. Do not induce vomiting. Call a physician immediately.

### 4.2. Most important symptoms and effects, both acute and delayed

Symptoms/effects after skin contact : Burns.  
Symptoms/effects after eye contact : Serious damage to eyes.  
Symptoms/effects after ingestion : Burns.

# ZirClean

## Safety Data Sheet

according to Regulation (EC) No. 1907/2006 (REACH) with its amendment Regulation (EU) 2020/878

### 4.3. Indication of any immediate medical attention and special treatment needed

Treat symptomatically.

## SECTION 5: Firefighting measures

### 5.1. Extinguishing media

Suitable extinguishing media : Water spray. Dry powder. Foam. Carbon dioxide.

### 5.2. Special hazards arising from the substance or mixture

Hazardous decomposition products in case of fire : Toxic fumes may be released.

### 5.3. Advice for firefighters

Protection during firefighting : Do not attempt to take action without suitable protective equipment. Self-contained breathing apparatus. Complete protective clothing.

## SECTION 6: Accidental release measures

### 6.1. Personal precautions, protective equipment and emergency procedures

#### 6.1.1. For non-emergency personnel

Emergency procedures : Ventilate spillage area. Avoid contact with skin and eyes. Do not breathe fume/mist/vapours.

#### 6.1.2. For emergency responders

Protective equipment : Do not attempt to take action without suitable protective equipment. For further information refer to section 8: "Exposure controls/personal protection".

### 6.2. Environmental precautions

Avoid release to the environment.

### 6.3. Methods and material for containment and cleaning up

Methods for cleaning up : Take up liquid spill into absorbent material.  
Other information : Dispose of materials or solid residues at an authorized site.

### 6.4. Reference to other sections

For further information refer to section 13.

## SECTION 7: Handling and storage

### 7.1. Precautions for safe handling

Precautions for safe handling : Ensure good ventilation of the work station. Avoid contact with skin and eyes. Do not breathe fume/mist/vapours. Wear personal protective equipment.  
Hygiene measures : Wash contaminated clothing before reuse. Do not eat, drink or smoke when using this product. Always wash hands after handling the product.

### 7.2. Conditions for safe storage, including any incompatibilities

Storage conditions : Store in a well-ventilated place. Keep cool.

### 7.3. Specific end use(s)

No additional information available

# ZirClean

## Safety Data Sheet

according to Regulation (EC) No. 1907/2006 (REACH) with its amendment Regulation (EU) 2020/878

### SECTION 8: Exposure controls/personal protection

#### 8.1. Control parameters

##### 8.1.1 National occupational exposure and biological limit values

No additional information available

##### 8.1.2. Recommended monitoring procedures

No additional information available

##### 8.1.3. Air contaminants formed

No additional information available

##### 8.1.4. DNEL and PNEC

No additional information available

##### 8.1.5. Control banding

No additional information available

#### 8.2. Exposure controls

##### 8.2.1. Appropriate engineering controls

###### Appropriate engineering controls:

Ensure good ventilation of the work station.

##### 8.2.2. Personal protection equipment

###### Personal protective equipment symbol(s):



###### 8.2.2.1. Eye and face protection

###### Eye protection:

Safety glasses

###### 8.2.2.2. Skin protection

###### Skin and body protection:

Wear suitable protective clothing

###### Hand protection:

Protective gloves

###### 8.2.2.3. Respiratory protection

###### Respiratory protection:

In case of insufficient ventilation, wear suitable respiratory equipment

###### 8.2.2.4. Thermal hazards

No additional information available

##### 8.2.3. Environmental exposure controls

###### Environmental exposure controls:

Avoid release to the environment.

### SECTION 9: Physical and chemical properties

#### 9.1. Information on basic physical and chemical properties

Physical state	: Liquid
Colour	: Deep Blue.
Appearance	: Gel.
Odour	: Odourless.

# ZirClean

## Safety Data Sheet

according to Regulation (EC) No. 1907/2006 (REACH) with its amendment Regulation (EU) 2020/878

Odour threshold	: Not available
Melting point	: Not applicable
Freezing point	: Not available
Boiling point	: Not available
Flammability	: Not applicable
Explosive limits	: Not available
Lower explosion limit	: Not available
Upper explosion limit	: Not available
Flash point	: Not available
Auto-ignition temperature	: Not available
Decomposition temperature	: Not available
pH	: > 13
Viscosity, kinematic	: Not available
Solubility	: Not available
Partition coefficient n-octanol/water (Log Kow)	: Not available
Vapour pressure	: Not available
Vapour pressure at 50°C	: Not available
Density	: Not available
Relative density	: Not available
Relative vapour density at 20°C	: Not available
Particle characteristics	: Not applicable

### 9.2. Other information

#### 9.2.1. Information with regard to physical hazard classes

No additional information available

#### 9.2.2. Other safety characteristics

No additional information available

## SECTION 10: Stability and reactivity

### 10.1. Reactivity

The product is non-reactive under normal conditions of use, storage and transport.

### 10.2. Chemical stability

Stable under normal conditions.

### 10.3. Possibility of hazardous reactions

No dangerous reactions known under normal conditions of use.

### 10.4. Conditions to avoid

None under recommended storage and handling conditions (see section 7).

### 10.5. Incompatible materials

No additional information available

### 10.6. Hazardous decomposition products

Under normal conditions of storage and use, hazardous decomposition products should not be produced.

## SECTION 11: Toxicological information

### 11.1. Information on hazard classes as defined in Regulation (EC) No 1272/2008

Acute toxicity (oral)	: Not classified
Acute toxicity (dermal)	: Not classified
Acute toxicity (inhalation)	: Not classified

# ZirClean

## Safety Data Sheet

according to Regulation (EC) No. 1907/2006 (REACH) with its amendment Regulation (EU) 2020/878

Potassium Hydroxide (1310-58-3)	
LD50 oral rat	333 mg/kg (Equivalent or similar to OECD 425, Rat, Male, Experimental value, Oral)

Skin corrosion/irritation : Causes severe skin burns.  
pH: > 13

Potassium Hydroxide (1310-58-3)	
pH	13.5 (0.56 %, 25 °C)

Serious eye damage/irritation : Assumed to cause serious eye damage  
pH: > 13

Potassium Hydroxide (1310-58-3)	
pH	13.5 (0.56 %, 25 °C)

Respiratory or skin sensitisation : Not classified  
Germ cell mutagenicity : Not classified  
Carcinogenicity : Not classified  
Reproductive toxicity : Not classified  
STOT-single exposure : Not classified  
STOT-repeated exposure : Not classified  
Aspiration hazard : Not classified

Potassium Hydroxide (1310-58-3)	
Viscosity, kinematic	Not applicable (solid)

### 11.2. Information on other hazards

No additional information available

## SECTION 12: Ecological information

### 12.1. Toxicity

Ecology - general : Before neutralisation, the product may represent a danger to aquatic organisms.  
Hazardous to the aquatic environment, short-term (acute) : Not classified  
Hazardous to the aquatic environment, long-term (chronic) : Not classified

Potassium Hydroxide (1310-58-3)	
LC50 - Fish [1]	80 mg/l (96 h, Gambusia affinis, Static system, Fresh water, Experimental value)

### 12.2. Persistence and degradability

Potassium Hydroxide (1310-58-3)	
Persistence and degradability	Biodegradability: not applicable.
Chemical oxygen demand (COD)	Not applicable
ThOD	Not applicable
BOD (% of ThOD)	Not applicable

### 12.3. Bioaccumulative potential

Potassium Hydroxide (1310-58-3)	
Bioaccumulative potential	Bioaccumulation: not applicable.

# ZirClean

## Safety Data Sheet

according to Regulation (EC) No. 1907/2006 (REACH) with its amendment Regulation (EU) 2020/878

### 12.4. Mobility in soil

#### Potassium Hydroxide (1310-58-3)

Surface tension	No data available in the literature
Ecology - soil	Low potential for adsorption in soil.

### 12.5. Results of PBT and vPvB assessment

No additional information available

### 12.6. Endocrine disrupting properties

No additional information available

### 12.7. Other adverse effects

No additional information available

## SECTION 13: Disposal considerations

### 13.1. Waste treatment methods

Waste treatment methods : Dispose of contents/container in accordance with licensed collector's sorting instructions.

## SECTION 14: Transport information

In accordance with ADR / IMDG / IATA / ADN / RID

### 14.1. UN number or ID number

UN-No. (ADR) : UN 1814  
UN-No. (IMDG) : Not applicable  
UN-No. (IATA) : Not applicable  
UN-No. (ADN) : Not applicable  
UN-No. (RID) : Not applicable

### 14.2. UN proper shipping name

Proper Shipping Name (ADR) : POTASSIUM HYDROXIDE SOLUTION  
Proper Shipping Name (IMDG) : Not applicable  
Proper Shipping Name (IATA) : Not applicable  
Proper Shipping Name (ADN) : Not applicable  
Proper Shipping Name (RID) : Not applicable  
Transport document description (ADR) : UN 1814 POTASSIUM HYDROXIDE SOLUTION, 8, III, (E)

### 14.3. Transport hazard class(es)

#### ADR

Transport hazard class(es) (ADR) : 8  
Danger labels (ADR) : 8



#### IMDG

Transport hazard class(es) (IMDG) : Not applicable

#### IATA

Transport hazard class(es) (IATA) : Not applicable

# ZirClean

## Safety Data Sheet

according to Regulation (EC) No. 1907/2006 (REACH) with its amendment Regulation (EU) 2020/878

### ADN

Transport hazard class(es) (ADN) : Not applicable

### RID

Transport hazard class(es) (RID) : Not applicable

#### 14.4. Packing group

Packing group (ADR) : III  
Packing group (IMDG) : Not applicable  
Packing group (IATA) : Not applicable  
Packing group (ADN) : Not applicable  
Packing group (RID) : Not applicable


#### 14.5. Environmental hazards

Dangerous for the environment : No  
Marine pollutant : No  
Other information : No supplementary information available

#### 14.6. Special precautions for user

##### Overland transport

Classification code (ADR) : C5  
Limited quantities (ADR) : 5I  
Excepted quantities (ADR) : E1  
Packing instructions (ADR) : P001, IBC03, LP01, R001  
Mixed packing provisions (ADR) : MP19  
Portable tank and bulk container instructions (ADR) : T4  
Portable tank and bulk container special provisions (ADR) : TP1  
Tank code (ADR) : L4BN  
Vehicle for tank carriage : AT  
Transport category (ADR) : 3  
Special provisions for carriage - Packages (ADR) : V12  
Hazard identification number (Kemler No.) : 80  
Orange plates :

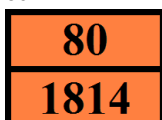

Tunnel restriction code (ADR) : E

##### Transport by sea

No data available

##### Air transport

No data available

##### Inland waterway transport

No data available

##### Rail transport

No data available

#### 14.7. Maritime transport in bulk according to IMO instruments

Not applicable



### SECTION 15: Regulatory information

#### 15.1. Safety, health and environmental regulations/legislation specific for the substance or mixture

##### 15.1.1. EU-Regulations

###### REACH Annex XVII (Restriction List)

Contains no substance(s) listed on REACH Annex XVII (Restriction Conditions)

###### REACH Annex XIV (Authorisation List)

Contains no substance(s) listed on REACH Annex XIV (Authorisation List)

###### REACH Candidate List (SVHC)

Contains no substance(s) listed on the REACH Candidate List

###### PIC Regulation (Prior Informed Consent)

Contains no substance(s) listed on the PIC list (Regulation EU 649/2012 concerning the export and import of hazardous chemicals)

###### POP Regulation (Persistent Organic Pollutants)

Contains no substance(s) listed on the POP list (Regulation EU 2019/1021 on persistent organic pollutants)

###### Ozone Regulation (1005/2009)

Contains no substance(s) listed on the Ozone Depletion list (Regulation EU 1005/2009 on substances that deplete the ozone layer)

###### Explosives Precursors Regulation (2019/1148)

Contains no substance(s) listed on the Explosives Precursors list (Regulation EU 2019/1148 on the marketing and use of explosives precursors)

###### Drug Precursors Regulation (273/2004)

Contains no substance(s) listed on the Drug Precursors list (Regulation EC 273/2004 on the manufacture and the placing on market of certain substances used in the illicit manufacture of narcotic drugs and psychotropic substances)

##### 15.1.2. National regulations

No additional information available

#### 15.2. Chemical safety assessment

No chemical safety assessment has been carried out

### SECTION 16: Other information

#### Indication of changes

Section	Changed item	Change	Comments
	Supersedes version of	Added	
	Revision date	Added	
	Issue date	Removed	
2.1	Classification according to Regulation (EC) No. 1272/2008 [CLP]	Modified	
2.2	Precautionary statements (CLP)	Modified	
7.1	Precautions for safe handling	Modified	
7.2	Storage conditions	Modified	

#### Full text of H- and EUH-statements:

Acute Tox. 4 (Oral)	Acute toxicity (oral), Category 4
Eye Irrit. 2	Serious eye damage/eye irritation, Category 2
H302	Harmful if swallowed.
H314	Causes severe skin burns and eye damage.

# ZirClean

## Safety Data Sheet

according to Regulation (EC) No. 1907/2006 (REACH) with its amendment Regulation (EU) 2020/878

### Full text of H- and EUH-statements:

H315	Causes skin irritation.
H319	Causes serious eye irritation.
Skin Corr. 1A	Skin corrosion/irritation, Category 1, Sub-Category 1A
Skin Corr. 1B	Skin corrosion/irritation, Category 1, Sub-Category 1B
Skin Irrit. 2	Skin corrosion/irritation, Category 2

Safety Data Sheet (SDS), EU

This information is based on our current knowledge and is intended to describe the product for the purposes of health, safety and environmental requirements only. It should not therefore be construed as guaranteeing any specific property of the product.